Hydrocolloid Dressing

Description

This Hydrocolloid Dressing contains carboxy-methylcellulose (CMC) which absorbs wound fluid and forms a gel to provide an optimal condition for wound healing. It is a top film that is bacteria-proof and waterproof.

Application of dressing

1. Wash hands before applying the dressing.
2. Clean and dry the wound and surrounding skin gently. Ensure no cream, ointment or similar product is on the skin.
3. Remove dressing from pouch.
4. Apply the dressing so that the sticky side sits directly over the wound.

When the product absorbs wound exudate, a white gel is formed and the appearance of the product will become white. The product should be changed when clinically indicated, when visible signs of exudate approach the edge of the product or after 7 days.

Removing the dressing:

5. Gently lift the corners of the product and remove it from the wound. Soaking the dressing with saline water where the wound is in contact can assist with removal.

Precautions for use

Seek medical advice immediately for all serious wounds and burns or if you experience redness or discomfort near the wound. A health care professional should frequently inspect and manage infected wounds, diabetic wounds and wounds which are solely or partially caused by arterial insufficiency, in accordance with local guidelines.

Dressing is sterile unless the individual pouch is opened or damaged. Do not use the dressing if the pouch is open or damaged. Single-use. Reuse of this single use product may create a potential risk to the user. Store in a cool, dry place (<35 degrees Celsius).

Intended Use

The intended purpose of our dressing is to provide an optimal condition for wound healing and create a barrier to protect wounds from bacterial and pathogens. Designed for dry to low exuding wounds. This is not intended to be a primary dressing used on chronic wounds for wounds involving the dermis.


Do not use

Re-use of the single use product may create a potential harm to the user. Reprocessing, washing, disinfection and/or sterilization may compromise product characteristics, causing additional risk of physical harm or infection to the user.

bpositive accepts no liability for any injury or loss that may arise if this product is used in a manner contrary to bpositive's current recommendations.

Basic UDI-DI 888501891021910BAV3

bpositive

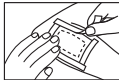Manufacturer:
BPositive Pte Ltd
www.bpositivehealthcare.com

Importer:
BPositive Pte Ltd
Singapore Business Federation Centre
160 Robinson Road, #14-04,
068914 Singapore

BPositive Healthcare Pty Ltd
Suite 18, 11 Wilson Street, South Yarra VIC,
3141 Australia

Local Authorised Representative:
T T MEDICAL MANAGEMENT SDN. BHD.
A-3A-03 EkoCheras, No:693, Batu 5 Jalan
Cheras, 56000 Kuala Lumpur, Malaysia
TL : +603 9134 5391

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BPOSITIVE PTE LTD,